



Hardware



Bridge Mount Quick Start Guide

Bridge Mount

Quick Start Guide for Bridge Mount for AI


Getting Started

What's in the box

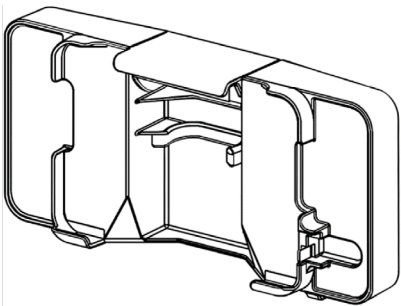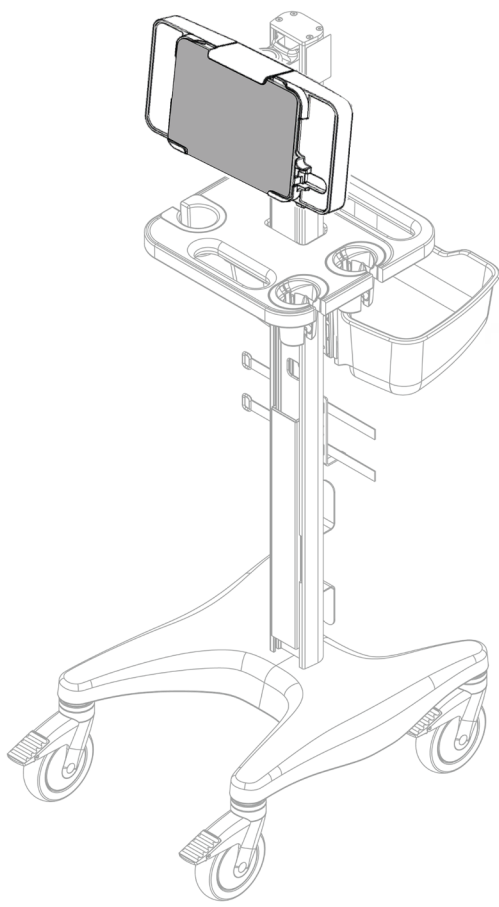

 (4) M4 x 8mm BHTS
(Button Head Torx Screw)









 (1) M3 x 8mm BHTS

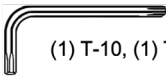 (1) T-10, (1) T-15, (1) T-20

 (2) Cable Clip


Installing tablet mount to column


Using a T20 Torx wrench, secure tablet mount with (4) M4 x 8mm BHTS as shown.

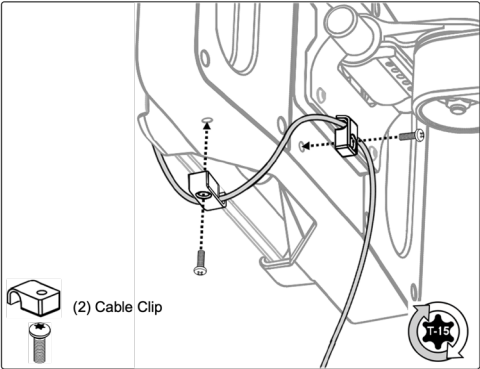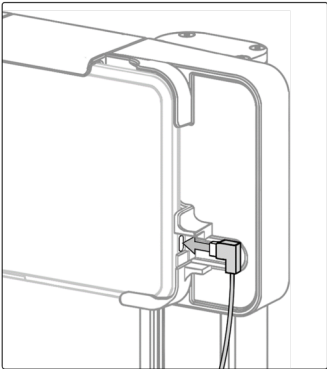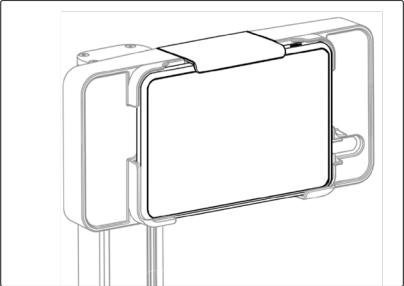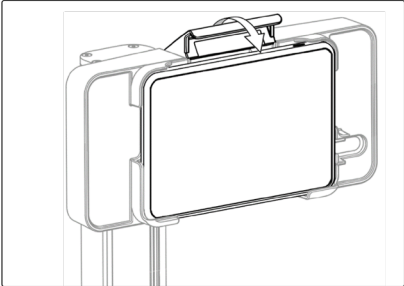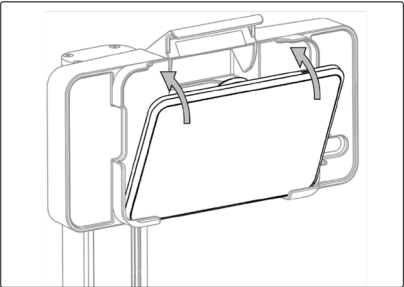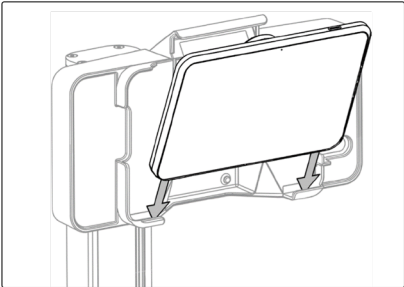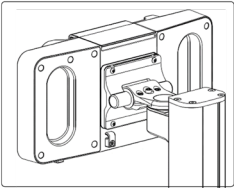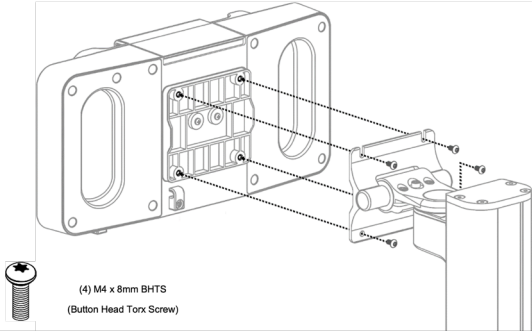

Bridge Mount

Quick Start Guide for Bridge Mount for AI


Installing Power Supply

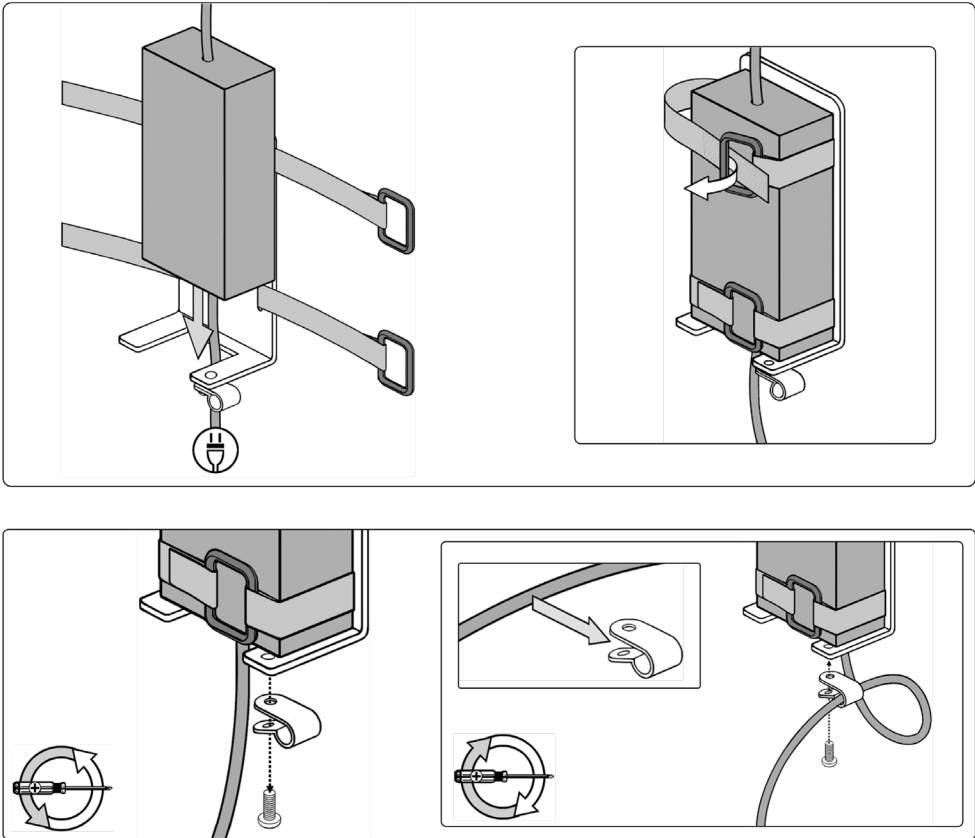

Cable routing

Column Cutouts

Route cable through cutouts in the column and up to tablet.





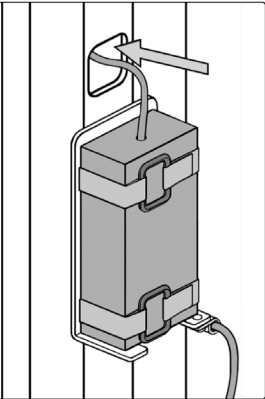


Cable Cleat

Store power supply cord on cable cleat.



Cable Clips

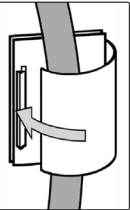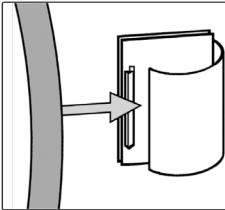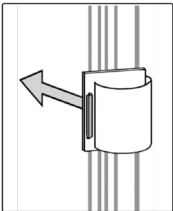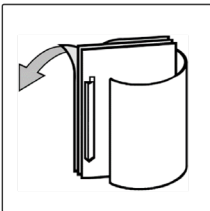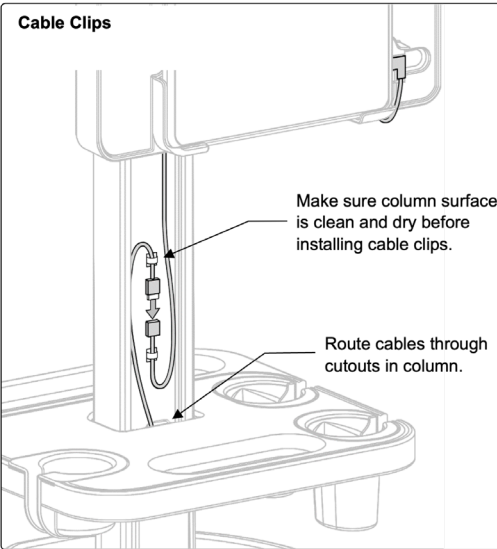



Bridge Mount

Quick Start Guide for Bridge Mount for AI

Tablet Release


To remove the Kosmos Bridge, press the release button (located on the back of mount). Pivot and lift the Kosmos Bridge out of mount. Power cable will automatically disconnect from Kosmos Bridge.

Release button is located on the backside within the right-side hand grip area.

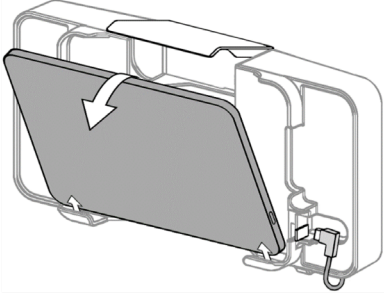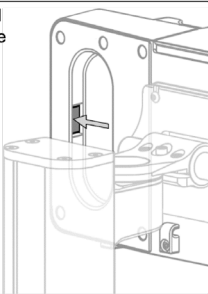

Locking the Kosmos Bridge

The Kosmos Bridge may need to be locked inside tablet holder. Once the Kosmos Bridge has been installed, use a T-10 Torx wrench to install (1) M3 x 8mm BHTS into the back of the tablet holder as shown below. This will lock the release button.

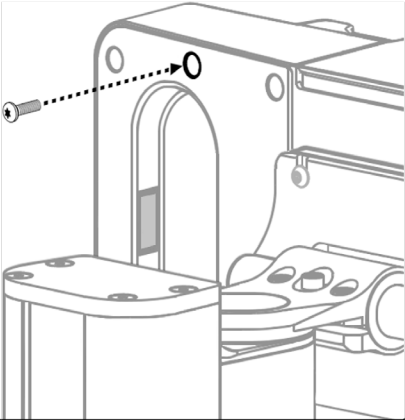

(1) M3 x 8mm BHTS





Customer Support



EchoNous, Inc.

8310 154th Ave. NE, Building B, Suite 200
Redmond, WA 98052

Call 1-844-854-0800
Email support@echonous.com



Android is a trademark of Google LLC.

© EchoNous, Inc., or its affiliates. All rights reserved.

EchoNous, echonous.com, and the EchoNous logo are trademarks of EchoNous, Inc. 8310 154th Ave. NE, Building B, Suite 200, Redmond, WA 98052.

Refer to user guide before use.

⚠ **Caution:** Federal (United States) law restricts this device to sale by or on the order of a physician.

Note: Not all features are available in all markets, please check with your local representative for availability in your region.

Please visit www.echonous.com for a current list of compatible tablets.

P007650-002_Rev A